



INSPIRE - 100

An Emergency Ventilator

TekMedika Pvt. Ltd.

(12) **United States Patent**
Nanda et al.








(10) Patent No.: US 12,465,716 B2
(45) Date of Patent: Nov. 11, 2025

INSPIRE-100 is a patent-pending ventilator design tailored for adult patients, providing comprehensive ventilatory support from initiation to weaning.

It features a user-friendly Human-Machine Interface with a clear front panel and advanced remote monitoring capabilities accessible via laptops, desktops, tablets, or mobile devices.

This device is particularly suited for facilities without access to compressed gas or oxygen pipelines and is built to withstand challenging environmental conditions.

INSPIRE-100 supports four commonly used ventilation modes: CMV, ACV, SIMV, and PSV. It offers a full range of breath parameters across all modes, ensuring complete assistance for patient-initiated breathing.

-  Unmatched Affordability
-  Unmatched Remote Monitoring via WiFi
-  Unmatched Ease-of-use
-  Works without Compressed Air Pipeline
-  Works with Oxygen Cylinder or Concentrator
-  Field Upgradeable
-  Complete Set of mainstream Parameters

Intended Use

Adult Patients only

- ✓ Acute Respiratory Distress Syndrome (ARDS)
- ✓ Chronic Obstructive Pulmonary Disease (COPD)
- ✓ Obese Hypoventilation Symptoms (OHS)
- ✓ Cheyne-Stokes breathing (CSR/CSA)
- ✓ Neuromuscular Diseases
- ✓ Pneumonia
- ✓ Asthma
- ✓ Drug Overdose
- ✓ Snake Bites

Technical Specs

Mode	Description
CMV	Continuous Mandatory
ACV	Synchronized Assist Control
SIMV	Synchronized Intermittent Mandatory
PSV (BiPAP)	Pressure Supported (BiPAP equivalent)

Volume Control Parameter	Range
Tidal Volume (ml)	200 - 600
Respiration Rate (bpm)	10 - 30
I:E Ratio	1:1 - 1:3
PEEP (cmH2O)	4 - 15
FiO2 (System Managed)	External

Pressure Support Parameter	Range
Support Pressure (cmH2O)	5 - 20
Flow Trigger Termination (%)	10 - 60


Full Set of ALARMS	
Max Pressure	Pressure Leak
Pressure Loss	Airway Blockage
Cough / Hiccups	System Temperature
Inconsistent Parameters	Extreme Parameters
Replace BVM	BVM Size
and many more ...	

Power Consumption 120W

Breath Synchronization for Patient Comfort


Breath Synchronization in SIMV Mode

e.g. RR=15 bpm (4 secs per breath)

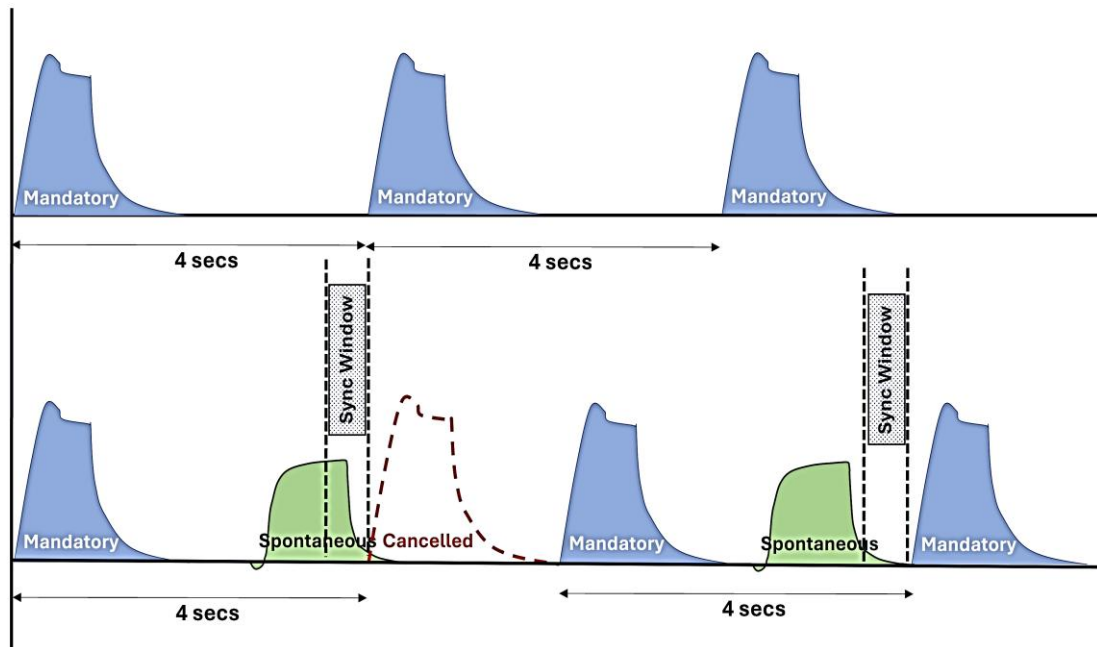


Breath Synchronization in ACV Mode

e.g. RR=15 bpm (4 secs per breath)


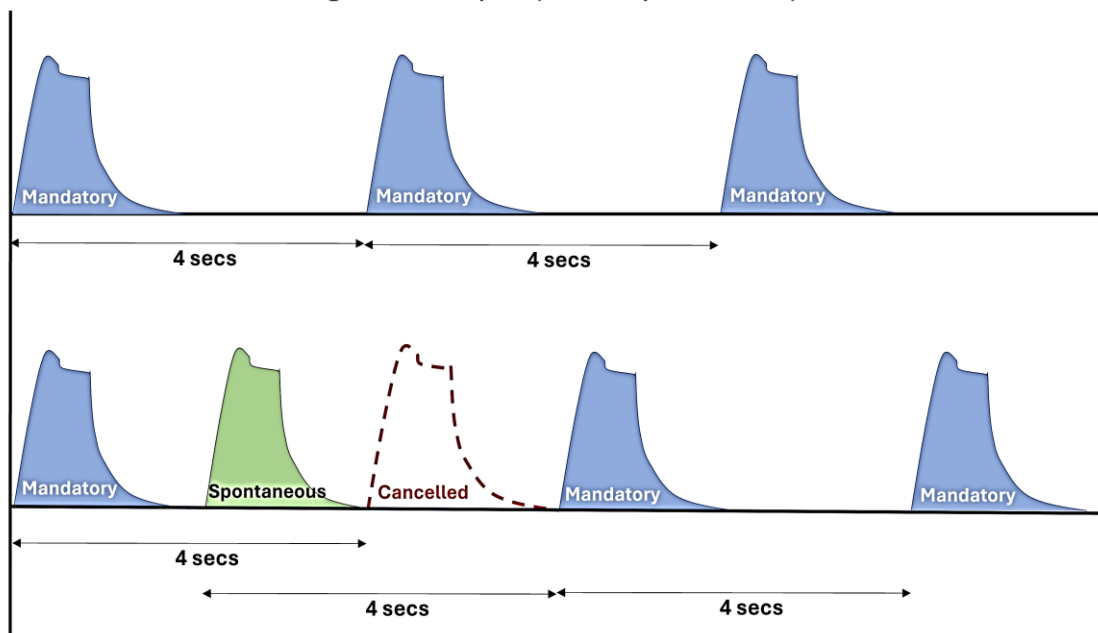


Breathing Circuit

Simple, Off-the shelf, Single-limbed and compatible with Standard Accessories

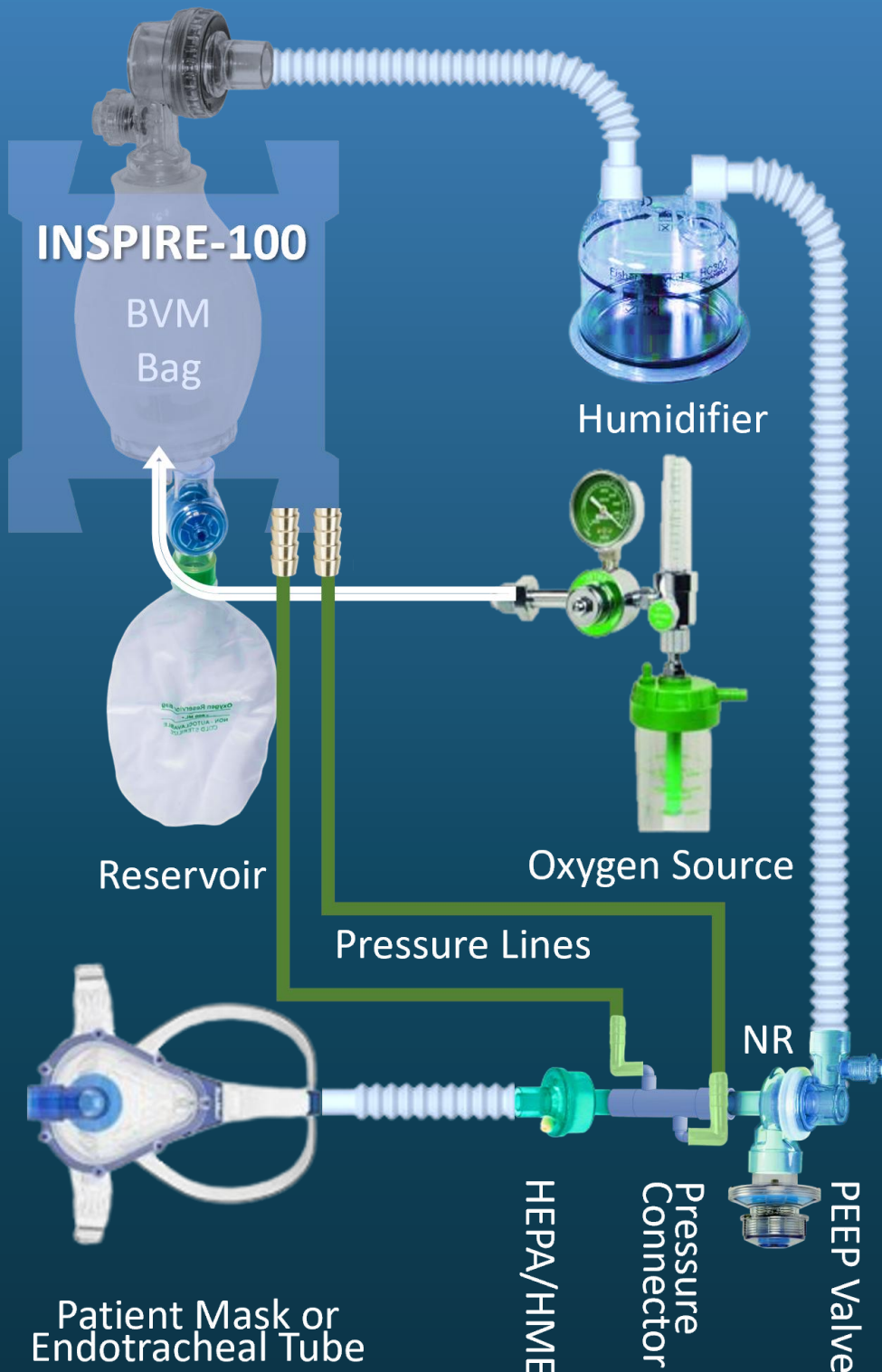

Remote Monitoring

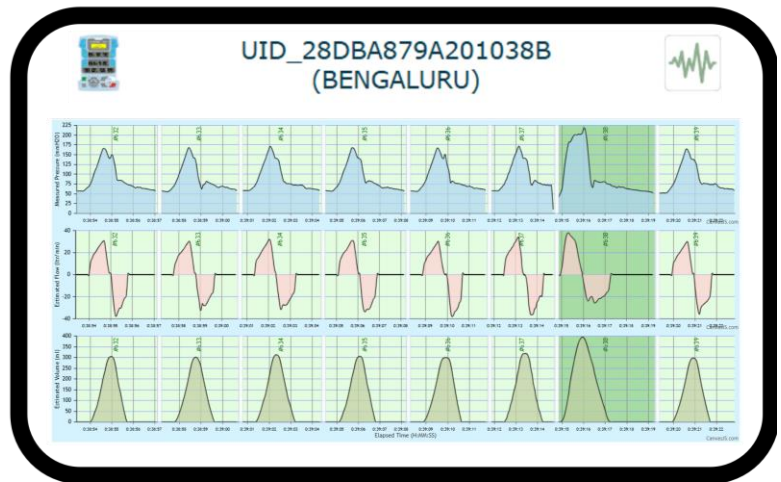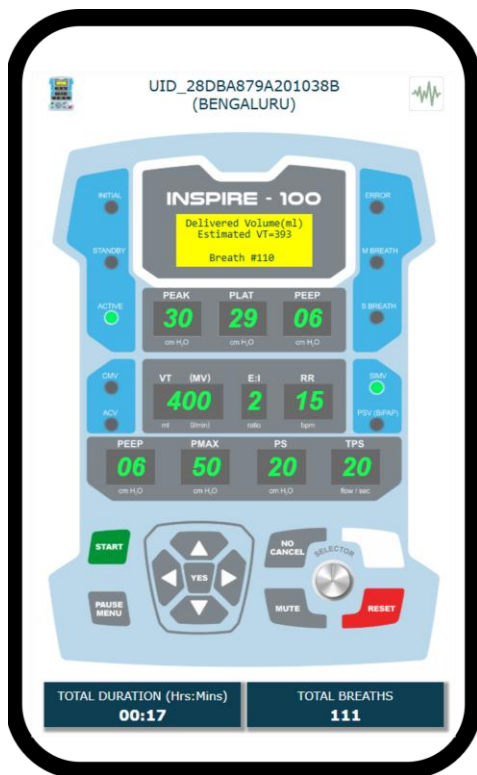

- ✓ Accessible on Laptops & Mobile devices
- ✓ Search & Range Selection
- ✓ Snapshots View
- ✓ Waveforms View
- ✓ Charts View
- ✓ Statistics View
- ✓ Recording and Playback
- ✓ Multi-system Display

Powerful Search for Combination of Events

1

128

Demo Long [09-05-2024 17:40:26]

MATCH CRITERIA

((MODE_SETTING == "SIMV") AND ((BREATH_TYPE == "MANDATORY") AND ((TIDAL_VOLUME < 250) OR (PEAK_PRESSURE > 30))))

MODE_SETTING

EQ

SIMV

AND

BREATH_TYPE

EQ

MANDATORY

AND

TIDAL_VOLUME

LT

250

ml

OR

PEAK_PRESSURE

GT

30

cmH2O

CLEAR

EXPORT

IMPORT

SEARCH

SEARCH RESULTS in Breath number range [1, 128]

Range	Matching Breaths			Parameter Values				
	Select	Number	Date	Time	MODE_SETTING	BREATH_TYPE	TIDAL_VOLUME	PEAK_PRESSURE
	<input type="checkbox"/>	9	09-05-2024	17:40:56	SIMV	MANDATORY	229	37
	<input type="checkbox"/>	10	09-05-2024	17:41:00	SIMV	MANDATORY	252	34
	<input type="checkbox"/>	11	09-05-2024	17:41:04	SIMV	MANDATORY	278	36
	<input type="checkbox"/>	43	09-05-2024	17:43:31	SIMV	MANDATORY	216	10
	<input type="checkbox"/>	44	09-05-2024	17:43:35	SIMV	MANDATORY	216	10
	<input type="checkbox"/>	46	09-05-2024	17:43:43	SIMV	MANDATORY	205	22
	<input type="checkbox"/>	47	09-05-2024	17:43:47	SIMV	MANDATORY	227	24
	<input type="checkbox"/>	48	09-05-2024	17:43:51	SIMV	MANDATORY	234	23
8 Matches: End of Search								


Range Selector to Navigate Any View Back and Forth in Time

1²³ Breath Range Selector


Selected Breath Range

	Breath Number	Breath Start Time		
		Day	Date	Time
From	520	Fri	26-Apr-2024	12:20:05
Upto	579	Fri	26-Apr-2024	12:22:58
SPAN	60	00:02:53		

Charts for All Parameters

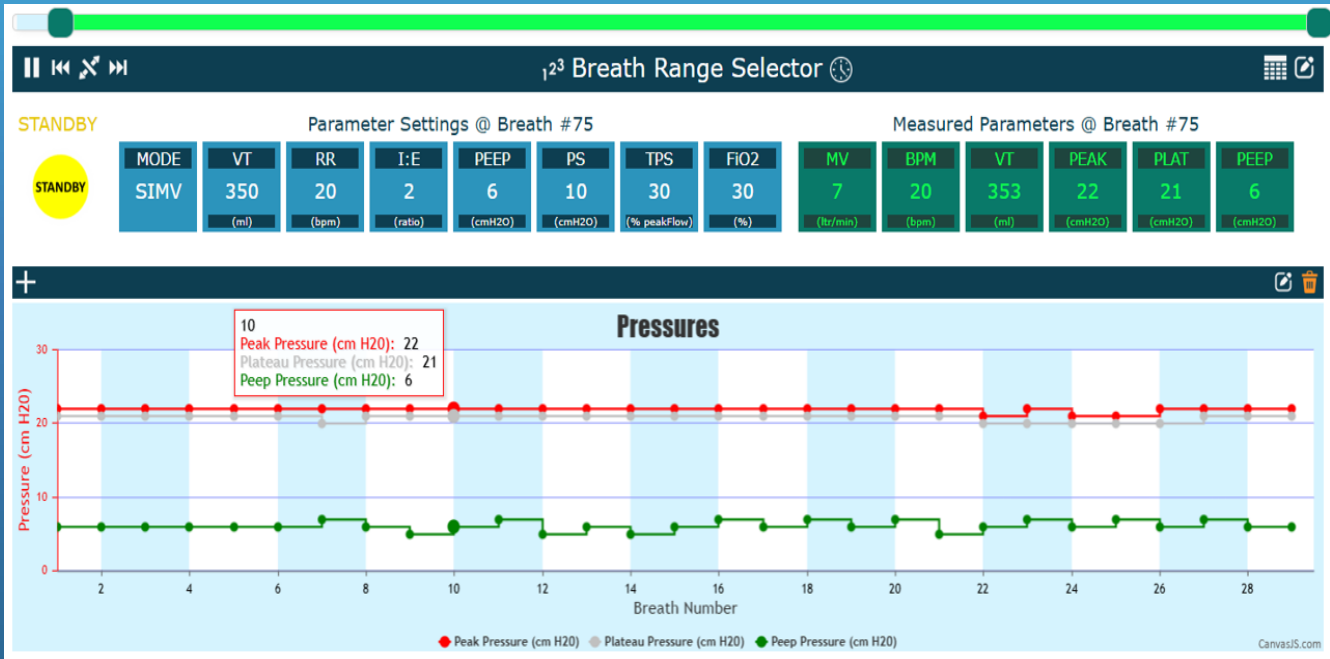


Pressure & Flow Waveforms

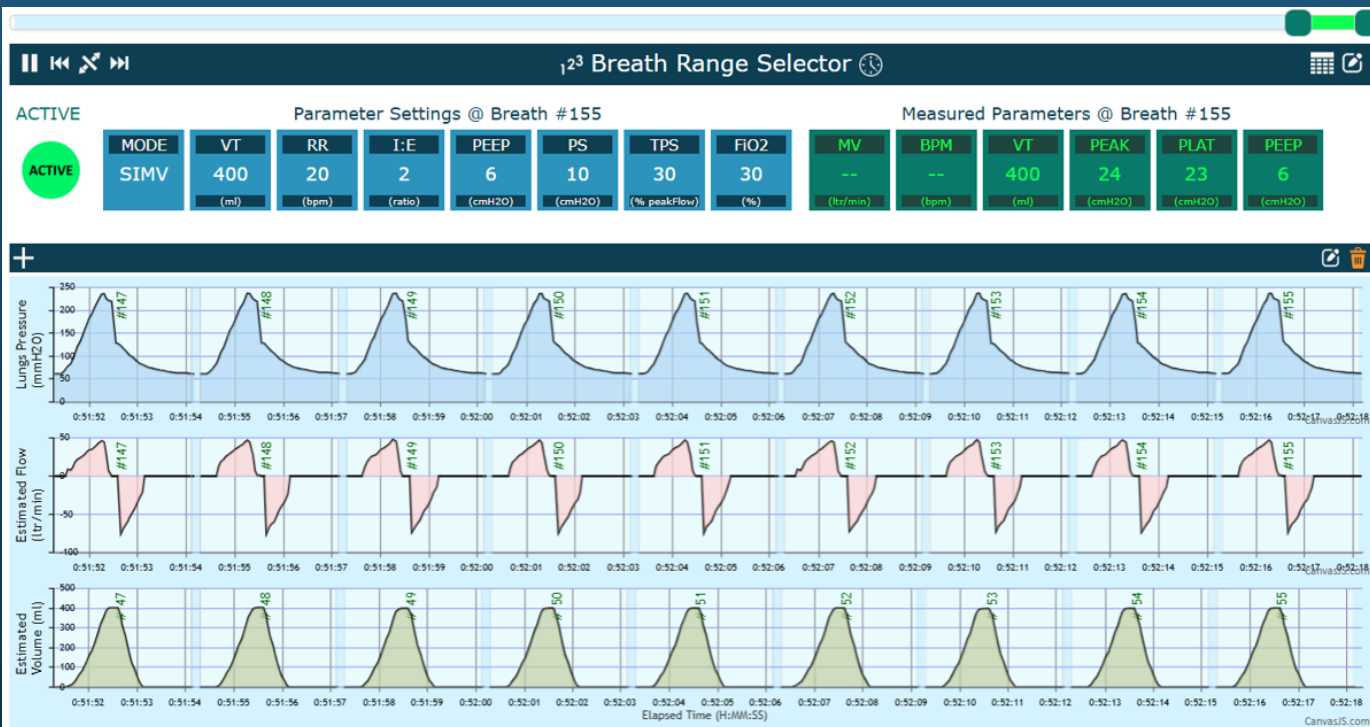

Normal Breaths			Maintenance Breaths			Exception Breaths	
Mandatory	Spontaneous	Spontaneous	Mandatory	PSV Mode	Spontaneous	Abnormal	Error
Volume Control	Volume Control	Pressure Control	Volume Control	Backup Breath	Volume Control	Seeming Breath	Causing Breath

Broken lines indicate few missing datapoints due to internet packet loss


Snapshots for each Breath

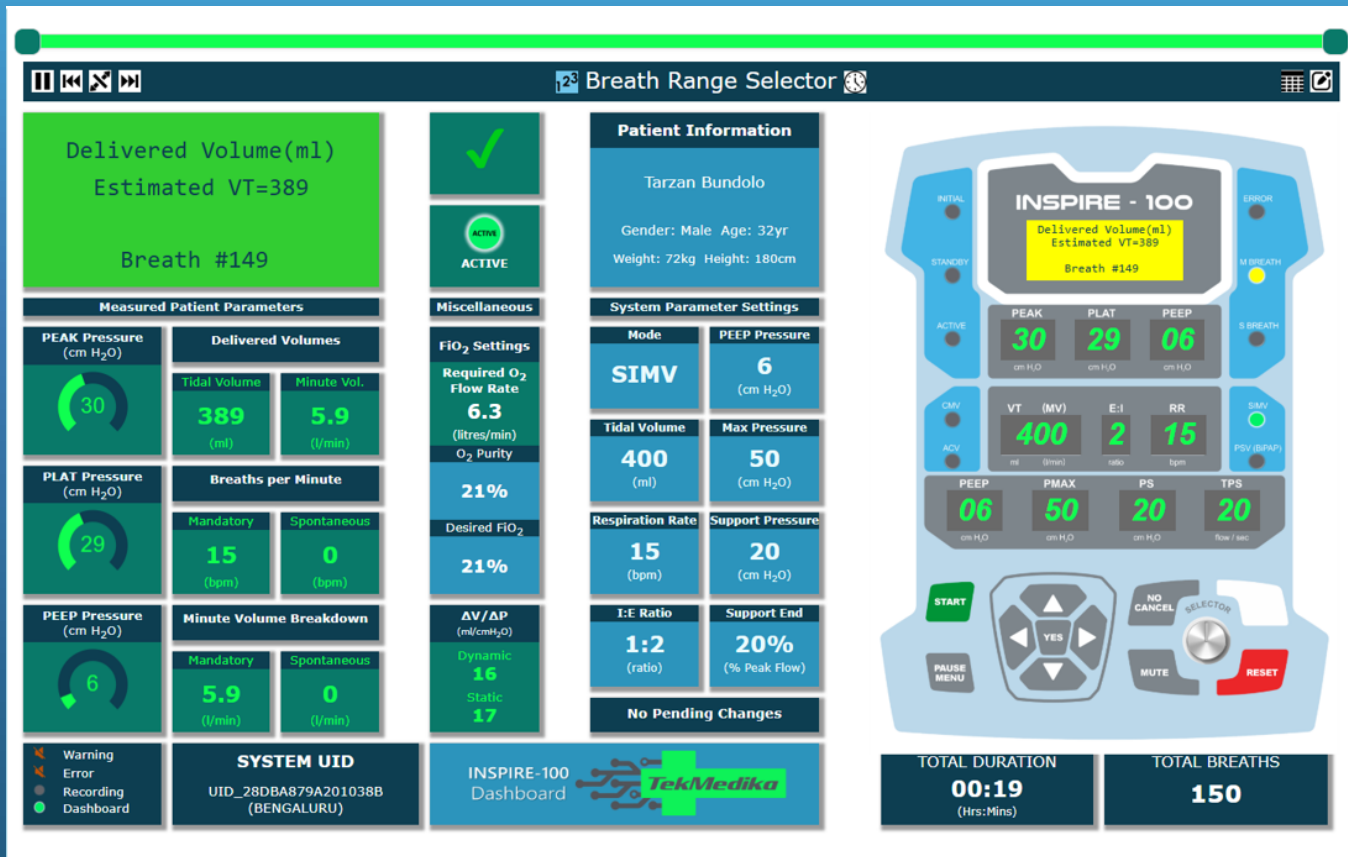

Detailed Statistics

Parameters Measured					Static Information																																			
Parameter	Units	Min	Max	Avg	<div>Patient Name: James Bond007 Gender: Male Age: 42yr Weight: 76kg Height: 185cm System Location: Bengaluru Location Altitude: 3000 ft Location Atmospheric Pressure: 930 cmH2O Location Atmospheric Oxygen: 19%</div>																																			
Peak Pressure	cmH2O	21	26	25.7																																				
Plateau Pressure	cmH2O	19	24	23.7																																				
PEEP Pressure	cmH2O	3	6	5.0																																				
Tidal Volume Delivered	ml	284	966	311.9																																				
Total Minute Volume	litres/min	0	5.8	4.3																																				
Mandatory Minute Volume	litres/min	0	4.6	4.1	<div>Parameter Settings Used</div> <table><tr><th>Parameter</th><th>Units</th><th>Values</th></tr><tr><td>Ventilation Mode</td><td>mode</td><td>SIMV</td></tr><tr><td>Tidal Volume</td><td>ml</td><td>300</td></tr><tr><td>Minute Volume</td><td>l/min</td><td>2</td></tr><tr><td>Respiration Rate</td><td>bpm</td><td>15</td></tr><tr><td>I:E Ratio</td><td>ratio</td><td>1:2</td></tr><tr><td>PEEP Pressure</td><td>cmH2O</td><td>5</td></tr><tr><td>Maximum Pressure</td><td>cmH2O</td><td>85</td></tr><tr><td>Support Pressure</td><td>cmH2O</td><td>20</td></tr><tr><td>Support Pressure Termination</td><td>%flow,secs</td><td>1.5</td></tr><tr><td>FIO2</td><td>%</td><td>35</td></tr></table>			Parameter	Units	Values	Ventilation Mode	mode	SIMV	Tidal Volume	ml	300	Minute Volume	l/min	2	Respiration Rate	bpm	15	I:E Ratio	ratio	1:2	PEEP Pressure	cmH2O	5	Maximum Pressure	cmH2O	85	Support Pressure	cmH2O	20	Support Pressure Termination	%flow,secs	1.5	FIO2	%	35
Parameter	Units	Values																																						
Ventilation Mode	mode	SIMV																																						
Tidal Volume	ml	300																																						
Minute Volume	l/min	2																																						
Respiration Rate	bpm	15																																						
I:E Ratio	ratio	1:2																																						
PEEP Pressure	cmH2O	5																																						
Maximum Pressure	cmH2O	85																																						
Support Pressure	cmH2O	20																																						
Support Pressure Termination	%flow,secs	1.5																																						
FIO2	%	35																																						
Spontaneous Minute Volume	litres/min	0	1.9	0.2																																				
Mandatory BPM	bpm	0	15	13.8																																				
Spontaneous BPM	bpm	0	2	0.2																																				
FIO2	%	35	35	35.0																																				
Static ΔV/ΔP	ml/cmH2O	0	18	15.9																																				
Dynamic ΔV/ΔP	ml/cmH2O	0	16	14.4																																				
System Temperature	degC	35	35	35.0																																				
Miscellaneous Information																																								
Information	Value																																							
Number of Breaths	120																																							
Number of Mandatory Breaths	118																																							
Number of Spontaneous Breaths	2																																							
Number of Maintenance Breaths	0																																							
Number of CMV Spontaneous Breaths	0																																							
Number of Missing Intervals (Packet loss)	0																																							
Number of WiFi Disconnects	0																																							
Number of Notifications	0																																							
Number of Warnings	0																																							
Number of Errors	0																																							
Sequence of Parameter Combinations																																								
MODE	VT/MV	RR	I:E	PEEP	PMAX	PS	TPS	FIO2	#BREATHS	Before#																														
SIMV	300	15	1:2	5	85	20	1.5	35	119	2																														

Multi-system Monitoring For Nurses Station

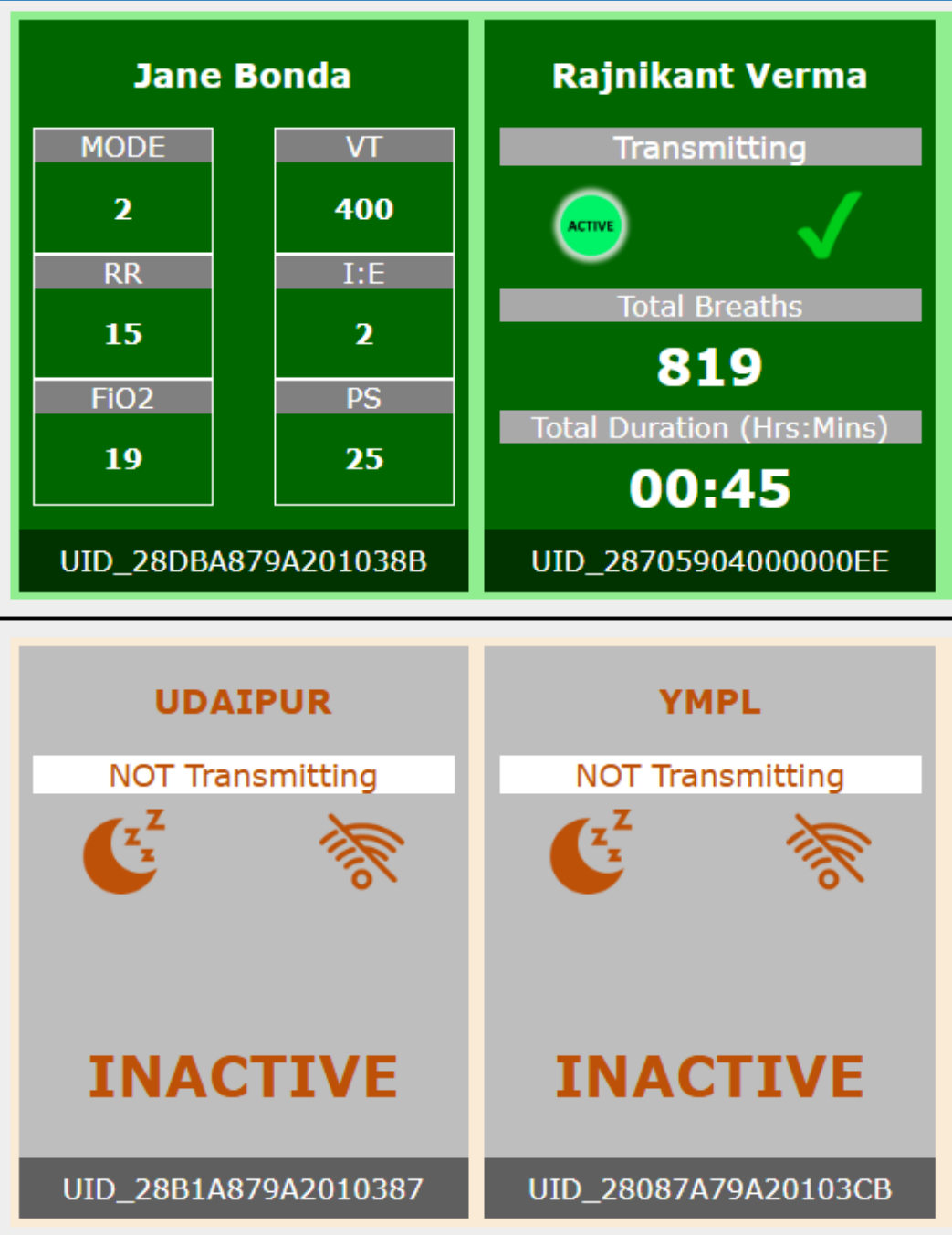

A summary state of all systems can be displayed on one screen. Each system is shown as a tile.



**Proudly
Designed
&
Made in India**